


Model 260-000 Bite Block System

Model 260-000 Bite Block System is a stand-alone patient head immobilizer system that needs no fixing to a treatment table and can be moved from room to room. The aluminum base plate has protruding pins which fit into the holes in the bottom of six different size head-neck supports, thus locating it to a fixed position every time it is being used. In some treatment setups it may be necessary to angle the head to the right or left. For these treatments, place the horizontal adjustable base on the base plate and lock in the position with two thumb screws. The head-neck supports also fit onto protruding pins located on the horizontal adjustable base plate or the pediatric adjustable base plate.

A three lock system allows for easy adjustment and reproducibility. Vertical and horizontal adjustment is accomplished using 15 cm scales and two locking knobs. The bite block tray angle is adjustable using a 180° scale and is locked with a knob.

An easy-to-remove individualized bite tray uses dental impression material. The impression procedure takes ten minutes using dental impression material (Model 266-100), available as an optional item. One pack will do 15 to 25 patients.

A fast setup can be accomplished by using the Nasal or Bridge trays. Rest the patient's head in the correct head-neck support, adjust the head for the desired chin to sternal stretch distance, lower the nasal or bridge tray, and lock.

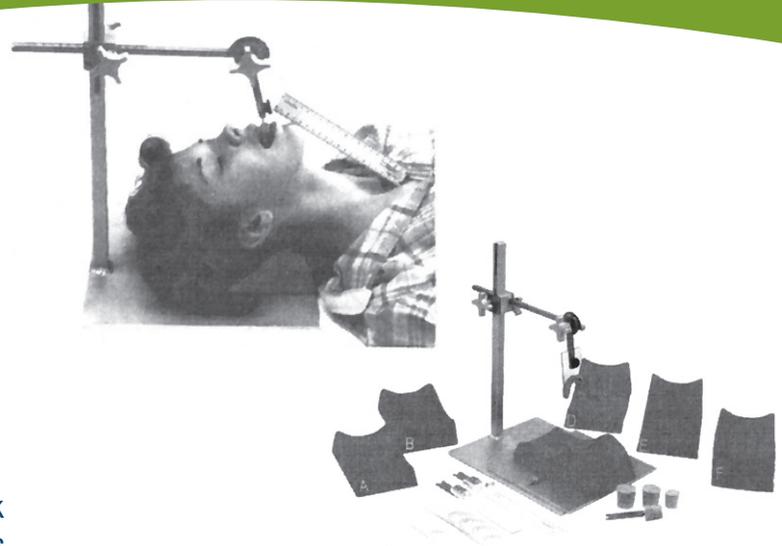

Features

- ▶ Keeps head in fixed reproducibility geometry
- ▶ Lightweight, portable, easy to use

Specifications

Size	43 cm vertical post welded to 25 cm x 33 cm base plate, aluminum base plate has an anti-skid rubber base
Lock system	Vertical, horizontal and tray angle
Scales	15 cm vertical, 15 cm horizontal and 180° tray angle
Shipping weight	7.6 kg (17 lbs)

Accessories

266-100	Dental Impression Material; one pack is enough for 15 to 25 patients
266-008	Pediatric Adjustable Base Plate

System includes:

1	Adjustable base plate
1	Bridge tray
6	Head-neck supports
1	Tongue tray
5	Adult bite trays
1	Chin tray
1	Child bite trays
1	20 cm white plastic rule
1	Nasal tray
1	Cork screw tray