


What's New from TeamBest?

Best® Total Solutions™ for Thermoplastics & Patient Immobilization

Huestis Medical CVO-2121 Thermoplastic Warming Oven

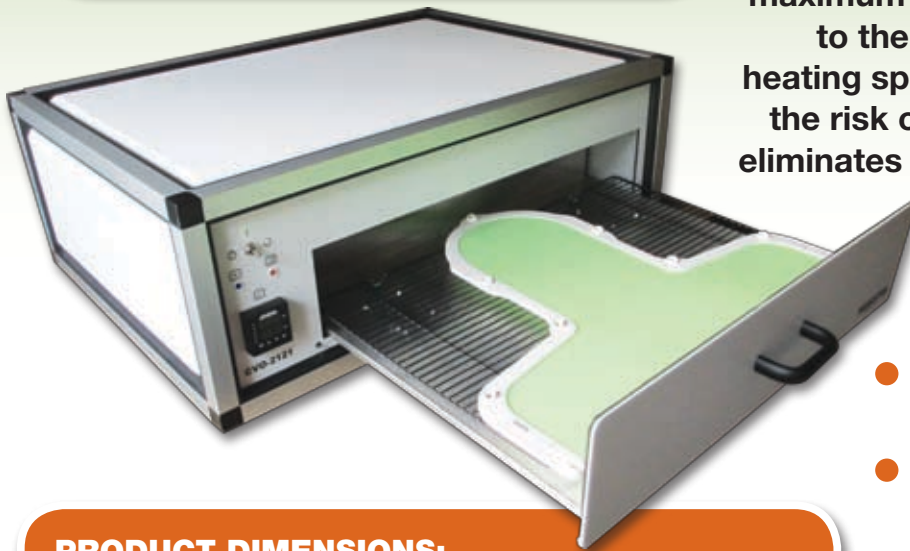

The Huestis Warming Oven softens thermoplastic material for radiation therapy to its optimum elasticity with safe dry heat, ensuring maximum comfort during application to the patient. Dry air convection heating speeds up workflow, reduces the risk of cross-contamination and eliminates other problems associated with water baths such as weight, spillage and humidity!

PRODUCT DIMENSIONS:

33.3" (85 cm) width x 25" (64 cm) depth* x 11.7" (30 cm) height
*The unit depth is 44" (112 cm) with drawer open

ELECTRICAL:

CVO-2121: 120 VAC, 14 Amps, Digital Control (NA market)
CVO-2221: 220 VAC, 7 Amps, Digital Control (EU market)



- Easy-sliding drawer opens fully to provide maximum access
- Multiple masks can be softened at the same time
- Can be used to mold head cushions without getting them wet
- Self-cleaning — a feature which reduces the potential for cross-contamination
- Can be used in CT rooms — no humidity that could otherwise effect the electronics of the imaging system
- Digital temperature setting in Celsius or Fahrenheit
- Non-stick surface
- Insulated case

© 2015 Best Medical International, Inc.

HUESTISMEDICAL
401.253.5500 huestismedical.com

CNMC+
615.391.3076 www.cnmc.co

Best® medical *international*
703.451.BEST (2378) teambest.com

Best® medical *canada*
613.591.2100 bestmedicalcanada.com

