


# Implant Needles



**A Bevel**


17g  
available in 10 cm, 15 cm, and 20 cm

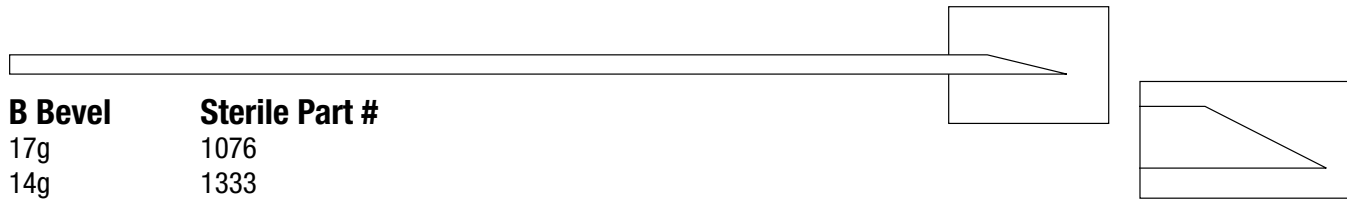

**B Bevel**

**Sterile Part #**

17g      1076  
14g      1333  
available in 10 cm, 15 cm, 20 cm, and 25 cm

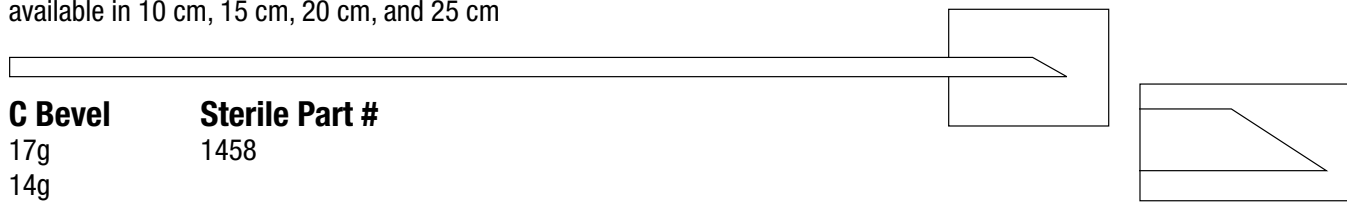

**C Bevel**

**Sterile Part #**

17g      1458  
14g  
available in 10 cm, 15 cm, 20 cm, and 25 cm  
\* 25 cm length only available in 17g

The Implant Needles are made of stainless steel and are hollow with different types of bevel at one end. Various lengths in different gauges are available, though the most commonly used needle is the 17g. Different types of Implant Needles are listed below:

	Part #	Part #	Part #
17g 10 cm	A—203—10	B—202—10	C—201—10
17g 15 cm	A—203—15	B—202—15	C—201—15
17g 20 cm	A—203—20	B—202—20	C—201—20
17g 25 cm		B—202—25	C—201—25
14g 10 cm		B—200—10	C—199—10
14g 15 cm		B—200—15	C—199—15
14g 20 cm		B—200—20	C—199—20

**Custom length is available upon request.**




**CAUTION:** This product should be carefully inspected by your staff prior to clinical use and should be used by authorized personnel only.

**Best Medical International, Inc.** 7643 Fullerton Road, Springfield, VA 22153 USA  
phone 703 451 2378 800 336 4970 fax 703 451 5228 www.bestmedical.com




# Template Needles


## Template Needles, 17g

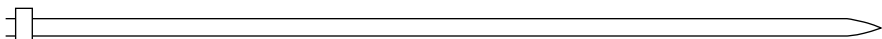With collar  Closed tip\*

Open Tip.....Part # 227  
Sterile Open End.....Part # 1395  
Sterile Closed End ...Part # 1397

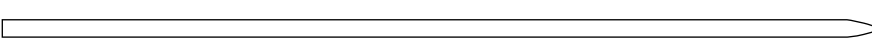Without collar  Open tip\*

.....Part # 196  
Sterile Open End.....Part # 1226  
Sterile Closed End ...Part # 1224

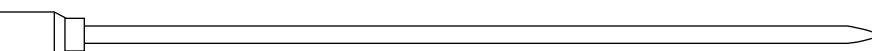With funnel  Open tip\*

Closed Sharp Tip .....Part # 231 Sterile.....Part # 1113  
Closed Tip.....Part # 198 Sterile.....Part # 1115  
Open Tip.....Part # 228

**\* All needles are available with open, closed or sharp tips  
Custom lengths available upon request; please specify length**

These 17 Gauge Needles are used with the templates, which improve the homogeneity of radiation. Needles of different lengths, closed end or open end, with or without collar, with or without funnel, with or without obturator are available. Catalog numbers of different template needles are given below:

OPEN END	STERILE	CLOSED END	STERILE
Tapered without Collar 10 cm, 15 cm, 20 cm, 25 cm... Part # 196.....Part # 1226		Sharp with Collar 10 cm, 15 cm, 20 cm, 25 cm... Part # 197.....Part # 1397	
Tapered with Collar 10 cm, 15 cm, 20 cm, 25 cm... Part # 227.....Part # 1395		Sharp without Collar 10 cm, 15 cm, 20 cm ..... Part # 193.....Part # 1224	
Template Needle w/Funnel & Open Tip 20 cm..... Part # 228.....Part # 1115		Template Needle (Pencil Point or Sharp) w/Funnel & Closed Tip 20 cm..... Part # 198.....Part # 1113	
Stylet..... Part # 332.....Part # 1228			



**CAUTION:** This product should be carefully inspected by your staff prior to clinical use and should be used by authorized personnel only.

**Best Medical International, Inc.** 7643 Fullerton Road, Springfield, VA 22153 USA  
phone 703 451 2378 800 336 4970 fax 703 451 5228 www.bestmedical.com




# Flexi-Needles

## Description

Flexi-Needles are sharp or blunt, closed-ended needle guides. They are made of biologically compatible material. They measure 20.2 cm from tip to end (custom lengths available upon request) and have nominal 0.3 cm long tip.

## Intended Use

Flexi-Needles are inserted transcutaneously, using special obturators. They are designed to be afterloaded with Iridium-192 or Iodine-125 ribbon sources. Two types of caps are available (Luer Cap or Slip Cap) to place over the loaded Flexi-Needle to hold the ribbons in place during irradiation. Simulator sources can be used for localization; obturators are left in place until the time each is afterloaded with the appropriate radioactive source.

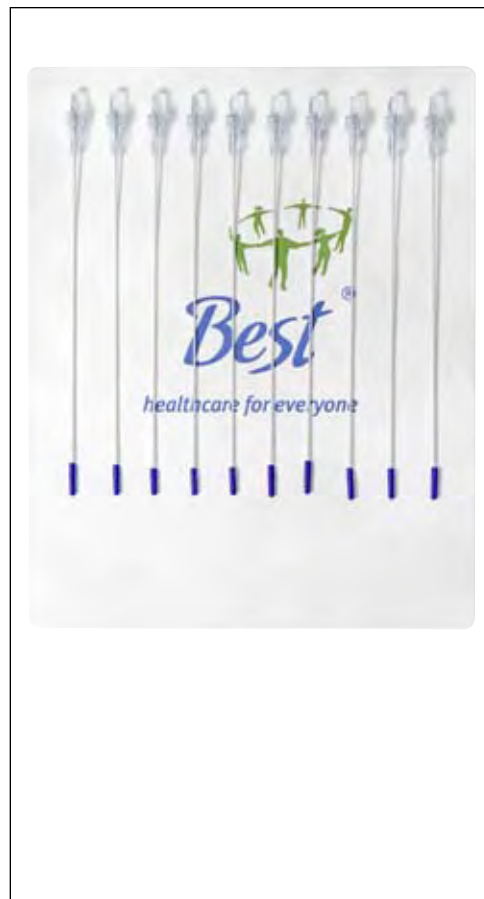

## WARNING

Obturators are made of Tungsten or stainless steel. Both obturator and Flexi-Needle may break if bone, hard skin, hard surfaces, or excessive forces are encountered.

## PRECAUTIONS

Caution: Device is shipped Non-Sterile or Sterile.

Caution: Flexi-Needles and Obturators are for single use only.

Caution: Federal law restricts this device to sale by or on the order of a physician.

Caution: Ensure Obturator is all the way to the tip of the Flexi-Needle prior to insertion.

Caution: Luer hubs may separate from the Flexi-Needle and Obturators if excessive force is applied.



**Best Medical International, Inc.** 7643 Fullerton Road, Springfield, VA 22153 USA  
phone 703 451 2378 800 336 4970 fax 703 451 5228 [www.bestmedical.com](http://www.bestmedical.com)



